

X-PRESS • PROFESSIONAL PRESSOTHERAPY

SELENIA ITALIA OUR HISTORY

SELENIA ITALIA SRL is a Beauty Company founded in 1991 by the Posarelli family, a leader in the field of professional aesthetics and aesthetic medicine. The company is specialized in the production and commercialization of innovative technologies and hitech cosmetics and is well known for its reliability, seriousness and maximum support to its customers. Selenia Italia is a certified company that distributes its products in Italy and in many European, African and Middle Eastern countries.



MADE IN ITALY

Innovation, effectiveness, reliability and Made in Italy design distinguish all of Selenia Italia's production, which has always been committed to the continuous pursuit of excellence. At the headquarters in Pisa, a specialized staff, in collaboration with a scientific team made up of doctors, engineers and technicians, devote themselves every day to the study and creation of new technologies and products to offer customers real, tangible and persistent results.

X-PRESS

EFFECTIVE DRAINAGE AND REMODELING OF THE LOWER LIMBS AND ABDOMEN

X-Press is an 8-way press therapy, ideal for counteracting the problems of body water retention and the imperfections of cellulite. It is indicated for drainage and microcirculation reactivation treatments, and for anti-cellulite and remodeling treatments.

Preset programs for every aesthetic problem, high quality bootleg with overlapping sectors, reliable electronics and many details make this device precious and indispensable.



X-PRESS


A PRECIOUS ALLY FOR YOUR BEAUTY CENTER!

It reproduces the caress and effects of an effective manual draining massage without the constant presence of the operator and allows to perform professional beauty treatments also combined with specific cosmetics and bandages in a practical and fast way: a great advantage for all beauticians!



8-SECTOR DE-LUXE BOOTLEG WITH STORAGE CLOTH

The bootleg, made of tear-resistant material, is fully sanitizable and adaptable in length and width to all body shapes through Velcro closures. The internal pressers can be replaced by the user without the need for technical maintenance. A practical container cloth to be placed on the bed contains both the bootleg and the connection tubes, allowing to neatly store the various components on the sides of the bed at the end of the treatment.

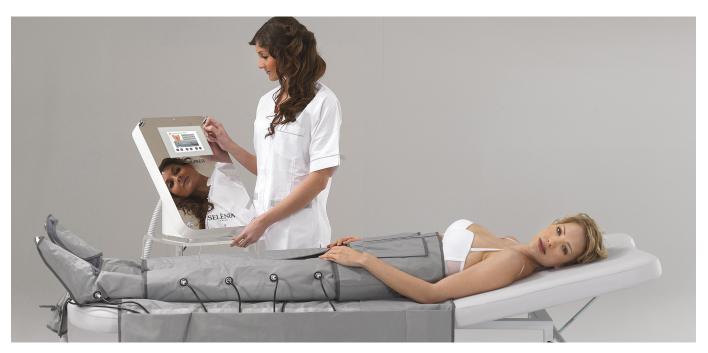





1 OPTIONAL

Jacket with 4 independent sectors for the treatment of the upper limbs.

2 OPTIONAL

Trolley for equipment.

MASSAGE WITH 4 DRAINAGE TECHNIQUES

X-PRESS allows to perform 4 sequential massage modes and various sequence maneuvers combined in 12 specific programs for the body and 4 programs for the upper limbs (optional). Each work program is designed to reproduce the manual massage best suited to the needs of the subject being treated. 4 drainage techniques allow to act on the body districts in a diversified and effective way.

DETOX SYSTEM

Specific drainage through the inflation of the sectors in correspondence with the main lymphatic ganglia: abdomen, groin, popliteal cavity, tibias.

POINT SYSTEM

Physiological drainage from bottom to the top. It consists in the sequential inflation of all the pressers, one sector at a time.

WAVE SYSTEM


It consists of alternating inflation, from bottom to top, of two sectors at the same time in order to increase the draining effectiveness without increasing pressure and without risks for blood capillaries.

LINE SYSTEM

It consists of inflating all sectors in sequence, from bottom to top, without releasing pressure until the end of the sequence, and then starting the cycle again. This innovative drainage system guarantees effective "compression" of the lower limbs, useful in cases of strong and deep retention.

WORKING SCREEN

An intuitive working screen allows the display of the SELECTED PROGRAM, the TREATMENT TIME, the PRESSURE and PRESSURE TIME parameters for each sector, and the DISPLAY of the ACTIVE SECTORS, instant by instant, throughout the session. The values set by the program can be changed at any time to be adapted to the customer's needs. Before each treatment, the device carries out a rapid PRE-INFLATION of all sectors to allow the bootleg and the pressers to adhere evenly to the whole body. At the end of the treatment, the EMPTYING button allows to quickly release the air contained in the pressers to facilitate the repositioning of the bootleg inside the container bag.

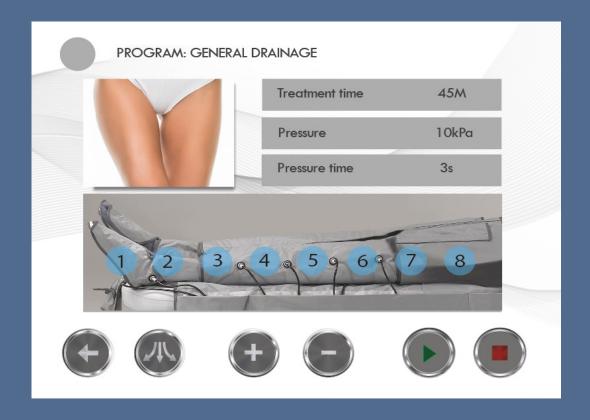

DURATION OF TREATMENT: from 30 to 50 minutes, depending on the area to be treated.

PROGRAMS

The pre-set programs allow the quick selection of the type of treatment to be carried out according to the customer's needs, combining 2 or more inflation techniques. Each program treats the body from the foot to the abdomen, but performs draining and stimulating maneuvers especially on the most critical areas and carrying out combinations of sequences designed to have the maximum draining effect in relation to the imperfection present.

GENERAL DRAINAGE

- PHYSIOLOGICAL
- LIGHT DRAINAGE
- INTENSE DRAINAGE
- STRONG DRAINAGE
- RELAXING

- REMODELING
- LIGHT LEGS
- REACTIVATING

CELLULITE

- COMPACT
- EDEMATOUS
- · SOFT
- CRITICAL POINTS

UPPER LIMBS

- LIGHT DRAINAGE
- INTENSE DRAINING 1
- INTENSE DRAINING 2
- INTENSE DRAINING 3



1.

DESIGN

Italian design, minimal, modern and elegant, that gives a touch of class to any beauty environment.

RESPECT FOR THE ENVIRONMENT

High quality electronic components and accessories. Equipment built with the minimum use of plastics to respect the environment. Eco-friendly case and trolley, made of recyclable material.



3.

5.7 " DISPLAY WITH HIGH RESOLUTION TOUCH SCREEN

A practical and intuitive interface with all the necessary information for customized programming and controls. Software available . in various languages.



BOOTLEG, PRESSERS, AND CONTAINER BAG

Bootleg in sturdy and washable tear-resistant material. 7 Sectors adaptable in length and width to all body shapes. Replaceable internal pressers. Practical container bag that allows to store all the accessories.

PRESET PROGRAMS

Massage programs for professional use of the device and specific for any type of body imperfection, even in combination with cosmetic products. 5.



INTEGRATED SELF-DIAGNOSTICS

A specific control system verifies the efficiency of the device in real time, reporting any anomalies or failures.

OUR RESULTS


BEFORE



AFTER



1 SESSION

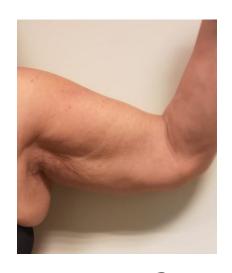



3 SESSIONS





SAFETY GUARANTEE

Device certified for Aesthetic use Compliant with inter-ministerial decree n ° 205 of 15/10/2015



SELENIA ITALIA QUALITY



QUALITY GUARANTEE

Selenia Italia is a company certified according to the EN ISO 13485-2016 Quality System..



Commercial name

General power supply

MAX absorbed power

Number of sectors

Internal protection

Category

Environmental operating characteristics

Humidity

Weight

Dimensions

X-PRESS

230V @ 50Hz

220 VA

8

Fuses n° 2 of 1,6AT (5x20mm)

Press-massage

Ambient temperature: 10 ° C to 40 ° C

30% to 75%

20Kg

35x52x33 cm

Classification: Device of type BF, class I

